


**OVERBED TABLE
SETUP GUIDE**

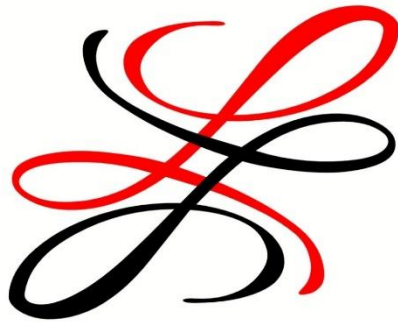

LJ CARE

1. OVERBED COMPONENTS

This Overbed Table comes with the following parts:

Description	Quantity
Tabletop with bracket mounted on bottom surface	x 1
Upright Post with quick adjustment lever and fixing knob	x 1
C-Shape Base Frame with 4 x black castor wheels (2 x locking)	x 1
14mm Bolts (Long) with washers and spacers (FOR BASE)	x 2
12mm Bolts (Short) with washers (FOR TABLETOP)	x 2

Thank you for choosing this OverbedTable.


This product has been designed to offer supreme stability, safety and adjustability.

It is important that you read this setup guide carefully before using your Overbed Table.

12 Month Warranty

LAMINATE TABLETOP

accommodates food trays, personal items, and provides a stable writing surface

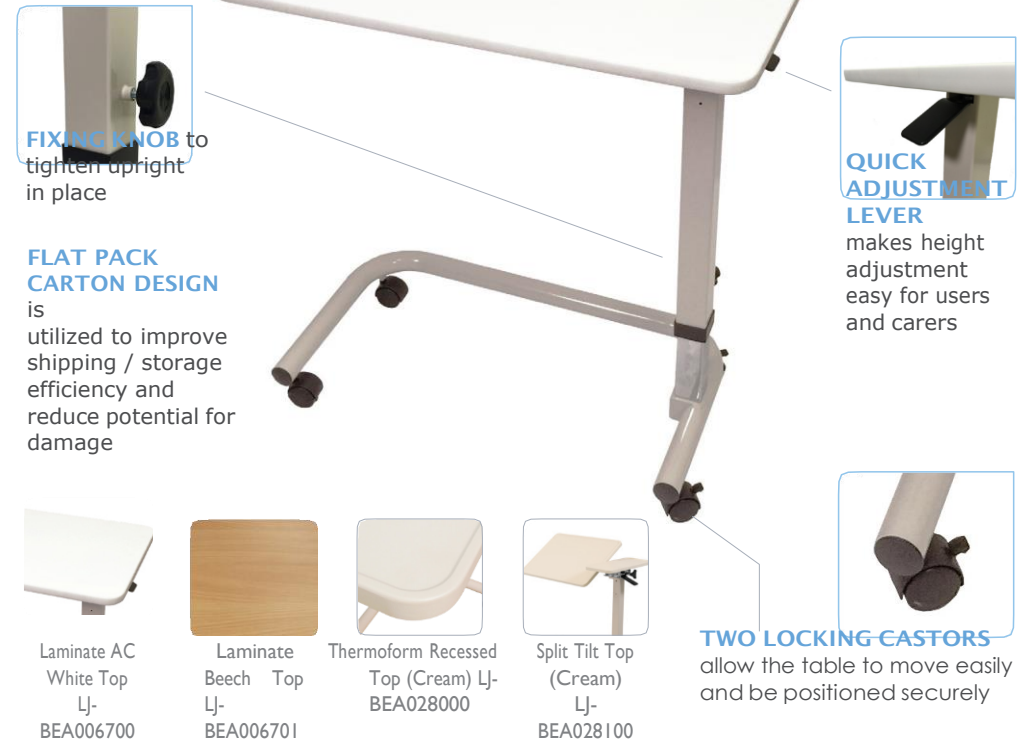

2. ASSEMBLY INSTRUCTIONS

Open the carton carefully. DO NOT USE a box cutter or sharp implement as this may damage the product. For OH&S, lay the carton on its side with the top of the carton open and slide out the Overbed Table contents.


1. Insert spacers, 14mm bolts and washers into the mounting plate on the C-shape base (from below).



2. Locate 14mm bolts into the base of the upright post and tighten. Ensure the upright post is straight on the base plate whilst tightening. Loctite or other fixing agents may be used if required.



3. Insert upright post table support bar into tabletop bracket and attach using 12mm short bolts and small washers.

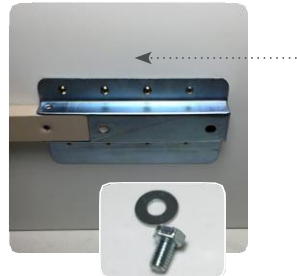

4. Remove the transportation screw from the upright and your over bed table is ready for use.



3. ADJUSTING YOUR TABLE

To Adjust your Over Bed Table:

1. Loosen the fixing knob on the upright post.



2. With one hand pull and release the quick adjustment lever and move the table up or down until the desired height is reached. Minimal force should be required.



3. Re-tighten fixing knob to lock in place.

Notes: